



CellinaPR
FRACTIONAL RF TECHNOLOGY

WORKSTATION

Primelaze®

We rise by our service

6 in One Workstation

This first-in-India workstation brings along both epidermal and dermal applications in one platform allowing professionals to work productively without changing any operational environment. Consisting of 6 different modalities, this platform has Microneedling RF, Subcision RF, RF Matrix, RF Thermolysis, Vaginal Rejuvenation, and RF Filler.



Cellina PR

India's only
Subcision RF
Device



01 Microneedling RF

Microneedles deliver RF energy directly to the targeted tissue in the skin and reform collagen & elastin effectively without any loss of energy unlike laser treatments.

Treatments

Acne Scars, Accidental, Traumatic & Suicidal Scars, Large Pores, Deep Wrinkles, Stretch Marks

Treatment Area: Middle & Deep Dermis




● Insulated Tip

● Non-Insulated Tip

Technology

User-friendly custom solution with advanced design


Insulated and non-insulated microneedles come in two types, allowing for tailored procedures without skin damage based on skin type and individual patients. Specially designed, incredibly thin gold-coated needles, along with upgraded motor, facilitate painless and fast procedures, as well as quick recovery in daily activities.

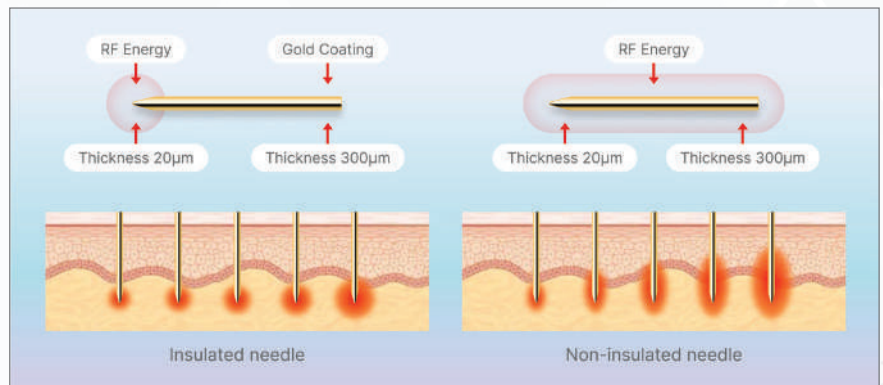

Precise and selective application

The dual wave combination of Cellina PR is capable of implementing various treatment modes to the maximum extent.

Continuous Wave: Irradiating the dermis layer with continuous energy generates high heat instantly, causing dermal layer coagulation and maximizing skin lifting and tightening effects.

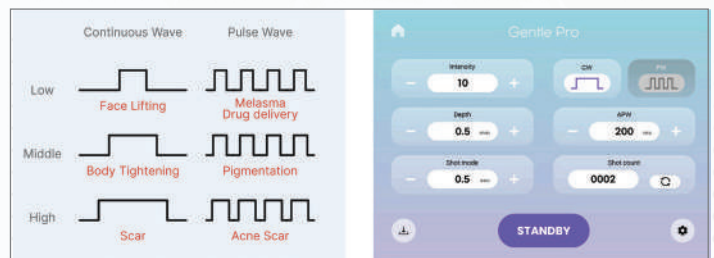

Pulse Wave: Various pulse modes deliver energy to the target area without damaging the surrounding tissue with minimal pain, effective for treatments such as pigmentation, acne, and melasma.

02 RF Matrix

The Fractional RF delivers thermal energy to the treatment area more safely, more effectively. It allows treatment regardless of skin type. It's safe for dark skin individuals, where use of lasers is avoided by practitioners.

Treatments

Skin Rejuvenation, Black Heads, Resurfacing, Fine Lines

Treatment Area: Epidermis & Superficial Dermis



03 Subcision RF


Treatments

Facial Contouring, Scars, Small Fat Reduction, Lifting, Facial Local Lypolysis



The principle is to help skin regeneration and to lead fibroblasts encouraged to produce collagen by itself through thermal energy disassemble irregular fiber cells under wrinkles & scars.

Patented cannula coating technology makes it easy for operator to deliver the RF energy at desired area. It generates RF energy only at the needle tip (1mm) to precisely deliver the energy in the tissue. FDA approved Parylene coating make the treatment safe.

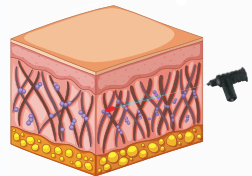

1. Approach to irregular layer of fibroblasts entangled under the wrinkles & Scars



2. The high-frequency heat energy delivered to the target to stimulate & decompose fibroblasts* Inducing collagen production & skin regeneration



3. Tangled fibroblasts are degenerated physically due to the high-frequency heat energy




4. Tangled fibroblasts are relocated to normal state due to the regenerative capacity of the skin & collagen production

04 RF Thermolysis

Treatments


Spider Vein, Acne, Open Pores, Syringoma, Papular Acne Scars

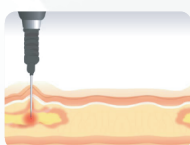


1. Deliver the high frequency heat energy to the target of acne




2. Decomposition & sterilization effect action to acne bacteria by thermal energy

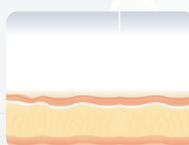


3. Remove acne through Cellina RF pen treatment

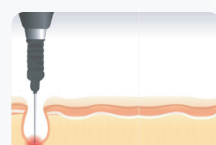

1. Deliver the high frequency heat energy to the pore

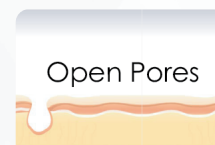

2. Pore shrinks by thermal energy

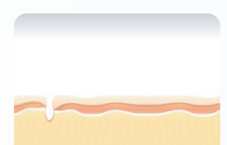

3. Pore shrink effect through the Cellina RF pen treatment



05 Vaginal Rejuvenation

Treatments

Vaginal Tightening, Stress Urinary Incontinence, Vulva Tightening

Three amazing changes with Safe & Simple 15 Min treatment.

Promotes proliferation of the collagen which makes the vaginal wall thickened & tightened

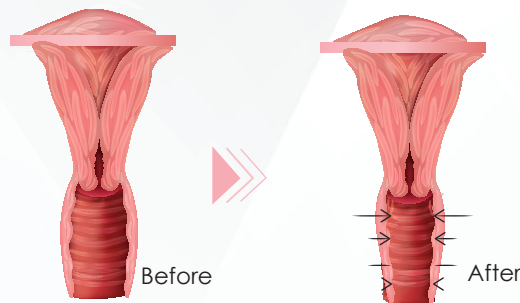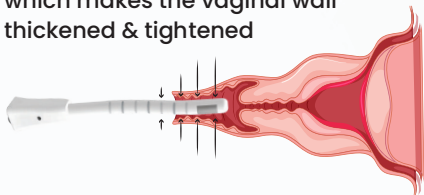

Cellina PR VT is the treatment for women. You can get good natural results for feminine rejuvenation from the 1st session. This treatment is different from other existing laser procedures, as it offers completely painless experience with no downtime. It offers thickening of vaginal uterine wall through superior collagen synthesis by Monopolar RF energy.





Vaginal Probe

06 RF Filler

Treatments

Wrinkle Treatment, Volume Augmentation, Facial Contouring

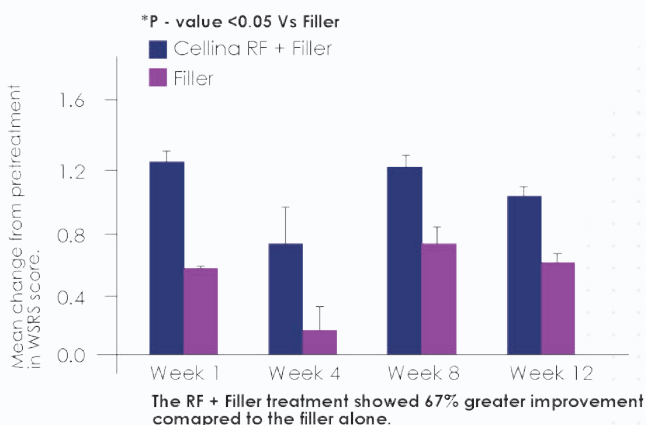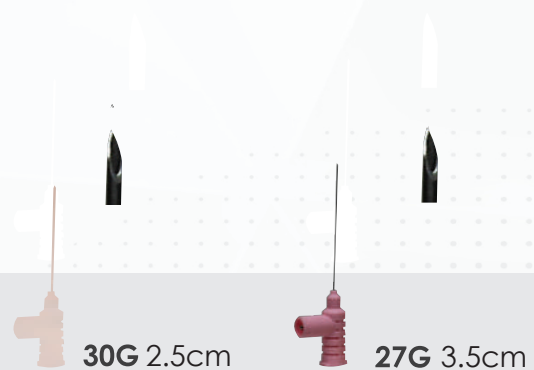

correction

RF tunneling improves filler effects & durability



Nasolabial folds based on 12 week observations



Specifications

Product Characteristics	Specifications
Technology	Radiofrequency
Frequency	1 MHz/2 MHz
Power	Max. 60W
Intensity	Up to 10 Levels
MNRF Needle Depth	Up to 4mm (0.1 mm Adjustable)
RF Exposure Time	100ms – 900ms
RF Cannula	23G, 21G, 19G & 15G
RF Delivery Modes	Continuous/ Pulse
Repetition Rate	0.5~3 Sec (0.5 Sec Step)
Display	Touch LCD Display
System Cooling	Air- Cooled
Weight	30Kg
Dimensions (Including Arm)	390mm(W) x 520mm(L) x 140mm(H)

Primelaze Private Limited

☎ 970 138 2034, 733 951 3001

✉ Email : Info@primelaze.com

🌐 www.primelaze.com

📍 Puducherry, Chandigarh, Mumbai

Manufactured by

Primelaze®

We rise by our service